



Shedding New Light On **MICROSCOPY**



Silicone Immersion Objective Series


Perfect for imaging 3D cell cultures and thick tissue samples. Silicone lenses enable visualization of cellular dynamics clearly and brightly when imaging at depth.

Research applications in neuroscience and cell biology using brain tissue, spheroids, organoids, and 3D cultures continue to push the limits of imaging in thick specimens. The need to image deep into such samples has never been more apparent. Nikon's silicone immersion objectives enable clear observation with high signal-to-noise deep into living tissue. The silicone immersion lens series features wide fields of view, high resolution, and evaporation-resistant oil, facilitating observations with ease.

With the addition of the new 60X objective, which employs newly developed glass for enhanced chromatic aberration correction, quantitative imaging in thick living specimens has never been more attainable.





Efficiently switch magnification for macro-to-micro imaging without changing the immersion medium

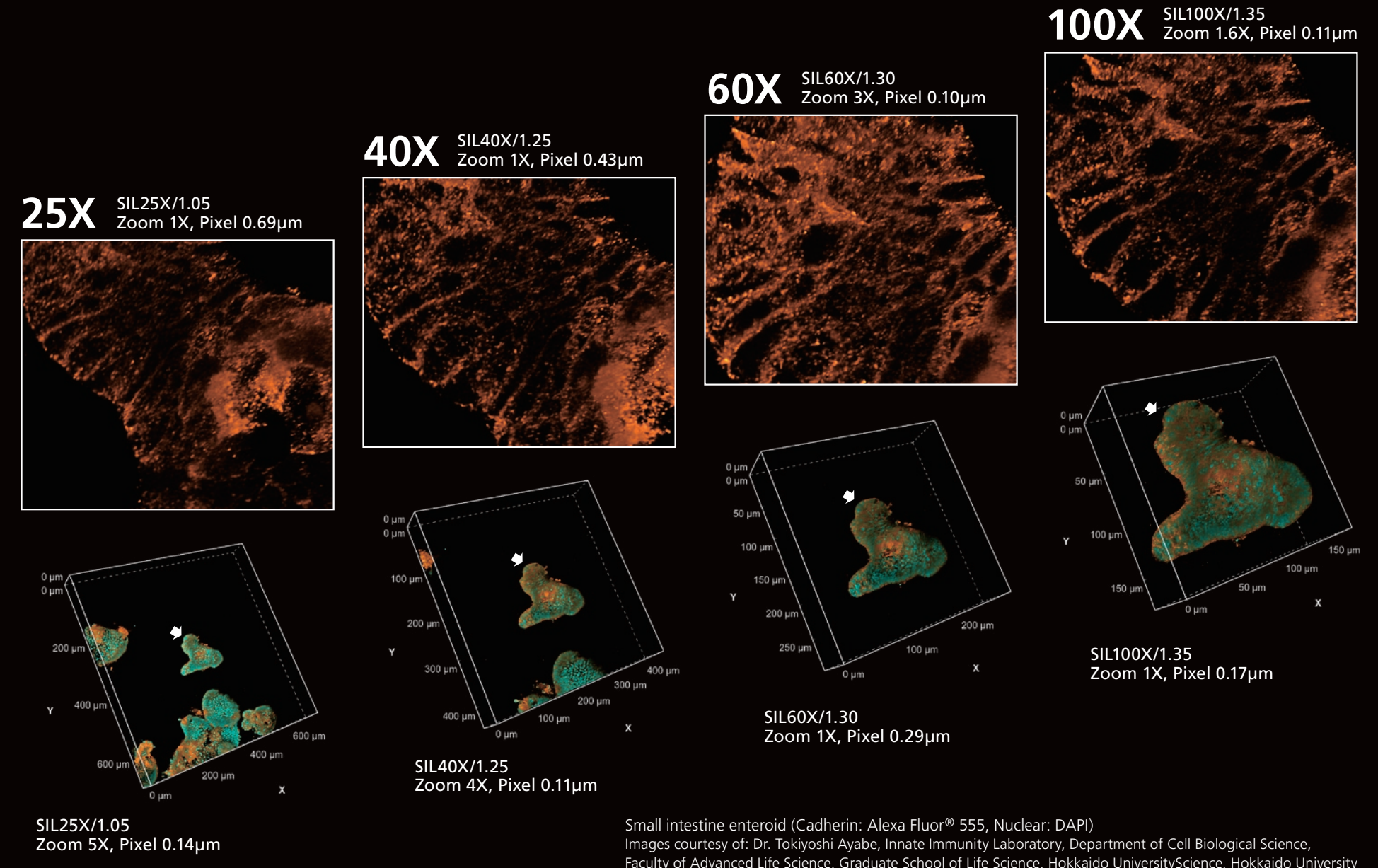


Small intestine enteroid (Cadherin: Alexa Fluor® 555, Nuclear: DAPI)
 Images courtesy of: Dr. Tokiyoshi Ayabe, Innate Immunity Laboratory, Department of Cell Biological Science, Faculty of Advanced Life Science, Graduate School of Life Science, Hokkaido University, Hokkaido University

Make high-resolution deep observations of thick specimens such as 3D cell cultures, organoids, and tissues

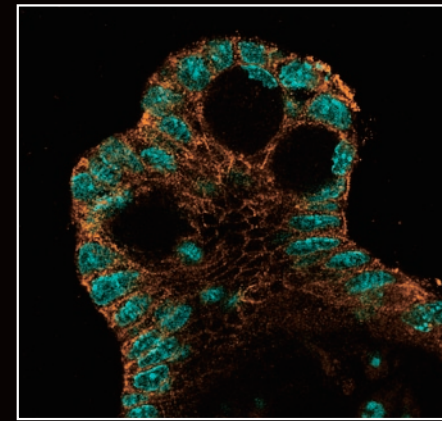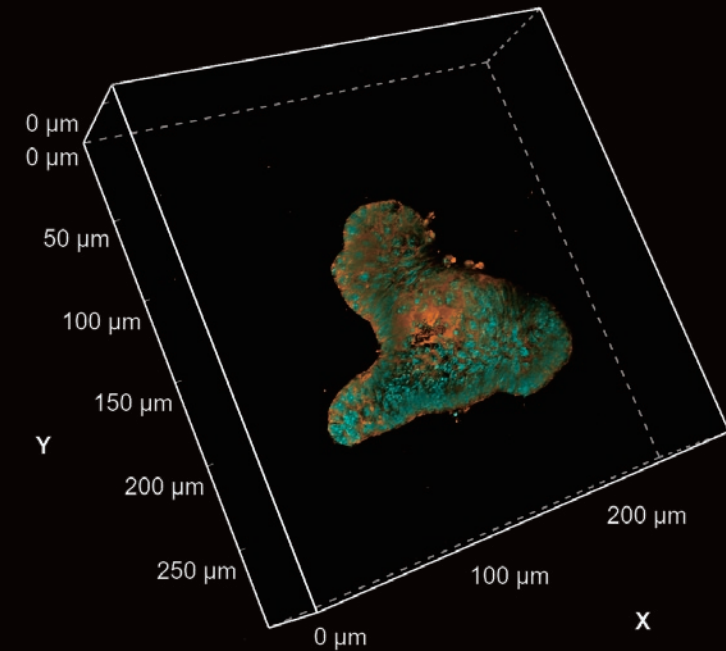

2D image surface

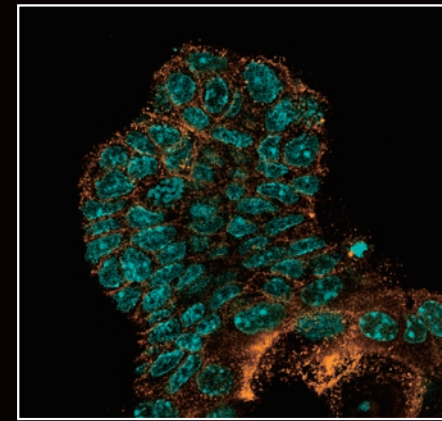


Deep 2D image

Small intestine enteroid (Cadherin: Alexa Fluor® 555, Nuclear: DAPI)
 Images courtesy of: Dr. Tokiyoshi Ayabe, Innate Immunity Laboratory, Department of Cell Biological Science, Faculty of Advanced Life Science, Graduate School of Life Science, Hokkaido University Science, Hokkaido University

Newly developed Short-wavelength Refractive (SR) glass


The new silicone immersion 60X objective employs high- and specialized-dispersion glass that was independently developed by Nikon and possesses extra low dispersion properties.

By refracting short-wavelength light to a higher degree, it is possible to collect a wider range of wavelengths, resulting in significantly enhanced chromatic aberration correction. In addition to axial chromatic aberrations, lateral chromatic aberrations can also be corrected.

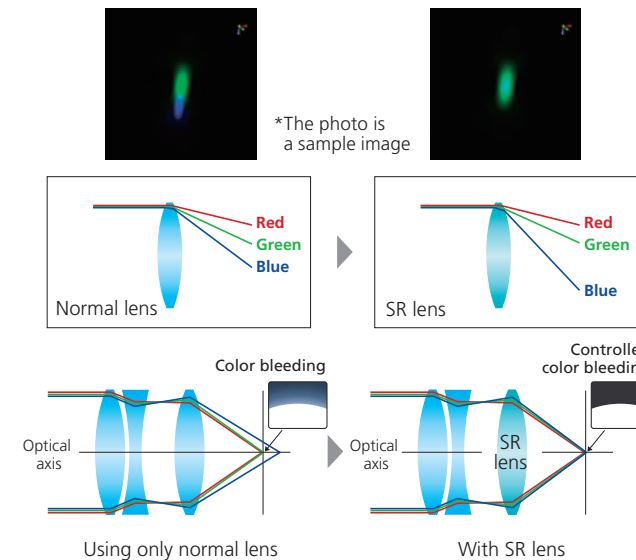


Observe live samples over long periods of time

It is also best suited for long time-lapse imaging of live cells without the immersion medium evaporating because silicone oil has low volatility even at 37°C.

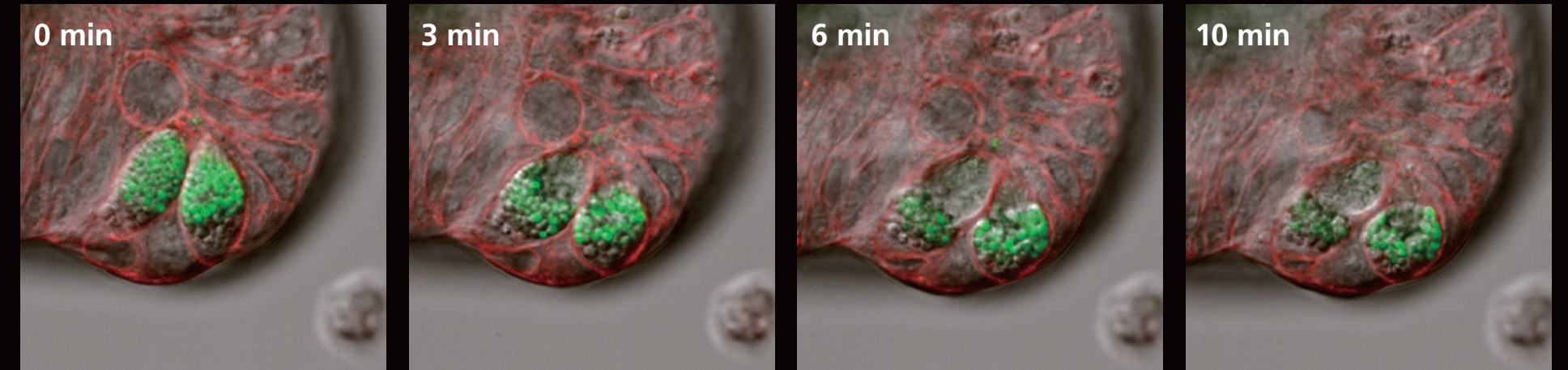

Time-lapse imaging of enteroid (25X) Images courtesy of: Dr. Yuki Yokoi, Dr. Kiminori Nakamura, and Dr. Tokiyoshi Ayabe, Innate Immunity Laboratory, Department of Cell Biological Science, Faculty of Advanced Life Science, Graduate School of Life Science, Hokkaido University

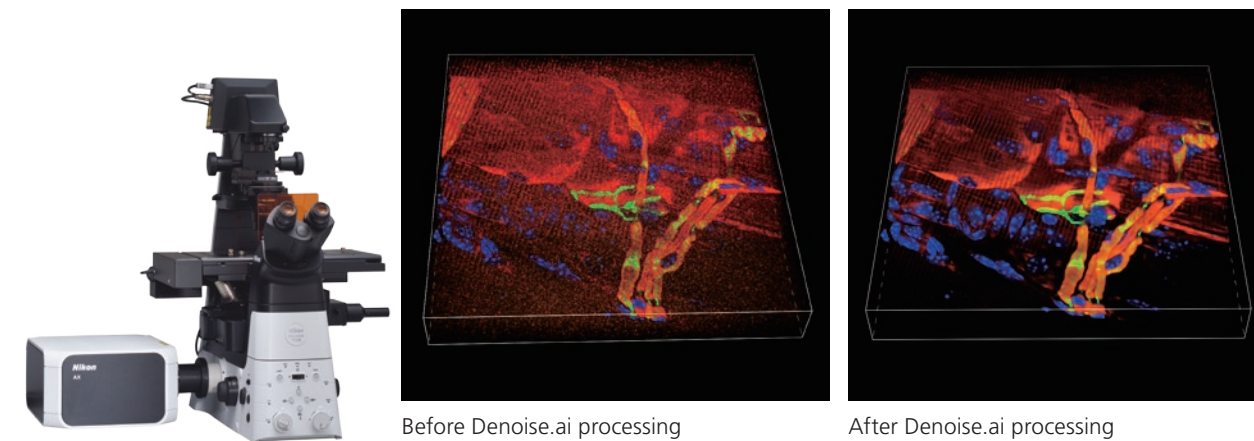

By using the Denoise.ai in conjunction with the AX confocal microscope, you can acquire even sharper images.

Details of examples are included in the application notes.

- High-speed, high-resolution 3D imaging using resonant scanner and Denoise.ai
- The Advantages of Resonant Scanning with Ultra Short Laser Exposure Times in Live Imaging
- In vivo Confocal Imaging of Mouse Organs that Clearly Captures Fast Dynamics (p-6)
- 3D Imaging of Intestinal Organoid (p-6)



Useful for a variety of research applications

Spheroid

Drug discovery

40X

Regenerative medicine

Small intestine organoid
Image courtesy of: Dr. Hidenori Akutsu and Dr. Tomoyuki Kawasaki of the Department of Reproductive Medicine, Center for Regenerative Medicine, National Center for Child Health and Development

Mouse nerves
Images courtesy of: Lin Daniel, PhD, SunJin Lab Co.

Fertilized mouse embryo
Image courtesy of: Dr. Yoshiteru Kai, Reproductive Medicine Research Center, YAMASHITA SHONAN YUME CLINIC

Growth

25X

Neuron

Drug discovery

60X

Neuron

100X


Dendrite

Neuron

Mouse synapse

Neutrophil flowing in blood vessel (time-lapse)
Images courtesy of: Professor Masaru Ishii, Department of Immunology and Cell Biology, Graduate School of Medicine, Osaka University

Comparison of immersion objective types

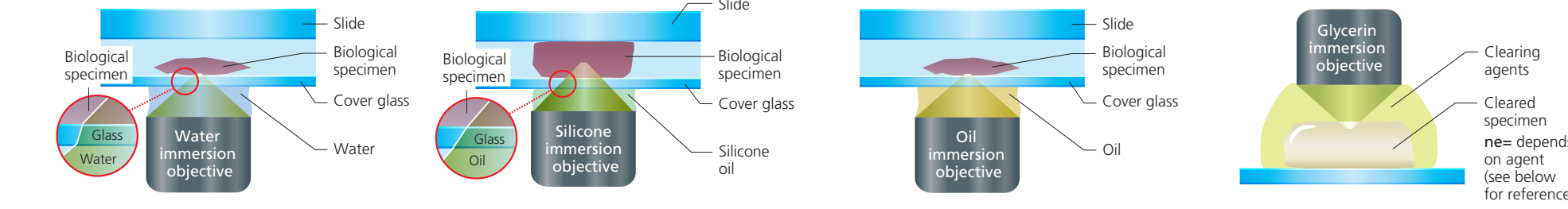

- Water immersion objective features**
 - Best for observing cultured cells and other specimens in culture solution because there is little spherical aberration as the refractive index is close to that of cultured cells and culture solution
 - Best for in vivo observations and live tissue fragments because it is highly biocompatible as the refractive index is close to that of biological fluids
 - Also effective with water immersion dispensers when it is necessary to supply immersion fluid, such as for observations of multiple points in multiple wells
- Silicone immersion objective features**
 - Little spherical aberration when making deep observations because the refractive index is close to that of many 3D cell cultures (spheroid, organoid, etc.) and organs with high refractive index (liver, etc.)
 - Best for long time-lapse observations because oil stable and resistant to drying even at 37°C.

Oil used for silicone immersion objective
Silicone immersion oil 30cc
- Oil immersion objective features**
 - Appropriate for observations that require high resolution because it has the highest NA of all immersion objective lenses
 - Best for observations of sample surfaces (vicinity of cover glass) because its refractive index is close to that of cover glass (ne=1.53)
 - Best for observations of anchored specimens that use a mounting medium because the refractive index is close to that of many mounting mediums
- Glycerin immersion (cleared specimen) objective lens features**
 - Possible to make deep observations of cleared specimens because glycerin has a refractive index close to many common clearing agents
 - Can use with a wide range of clearing agents
 - Can optimize spherical aberration correction to match the refractive index of the clearing agent using the correction collar

Examples of clearing agents appropriate for glycerin immersion (for clearing) objective lens
Sca/e(n=1.38)*, CUBIC2(n=1.48), LUCID(n=1.47), FocusClear(n=1.45), CRARITY(n=1.45), SeeDB(n=1.48)

Comparison of immersion objective features

Immersion medium (possible magnification range)	Water (20X–60X)	Silicone (25X–100X)	Oil (40X–100X)	Glycerin (10X–20X)
Immersion liquid refractive index	Water (ne=1.33)	Silicone (ne=1.41)	Oil (ne=1.52)	Clearing agent (nd=1.33–15.1; 10X, nd=1.44–1.50; 20X) Glycerin (ne=1.45)
Appropriate live samples for observation	Cultured cells Animal tissue and organs (mouse brain, etc.) with low refractive index Observations that must be biocompatible (in vivo observations, etc.) Observations that require multi-point observations of multiple wells (use of water dispensers)	3D cell cultures (spheroid, organoid, etc.) Animal tissue and organs with high refractive index (liver, etc.) Long time-lapse observations	Cultured cells and tissue fragments that require high resolution (sample surface)	—
Appropriate anchored samples for observation	Anchored cultured cells Anchored organ fragments with low refractive index (mouse brain, etc.)	Anchored fragments, such as organs with high refractive index (liver, etc.)	Cultured cells (with mounting medium) Tissue fragment (with mounting medium)	Cleared samples (spheroid, organoid, organs, etc.)

* Refractive index depends on the organism

Silicone immersion objective specifications

Model	Immersion	NA	W.D. (mm)	Cover glass thickness	Correction ring	Bright field	Dark field	DIC	Phase contrast	Polarizing	Fluorescence		Ti2-E PFS
											UV	Visible light	
CFI Plan Apochromat Lambda S 25XC Sil	Sil	1.05	0.55	0.11–0.23	✓	✓	Oil	✓		✓	✓	✓	
CFI Plan Apochromat Lambda S 40XC Sil	Sil	1.25	0.30	0.13–0.21 (23°C) 0.15–0.23 (30°C)	✓	✓	Oil	✓		✓	✓	✓	
CFI Plan Apochromat Lambda S 60XC Sil	Sil	1.30	0.30	0.15–0.19	✓	✓		✓		✓	✓	✓	
CFI SR HP Plan Apochromat Lambda S 100X Sil	Sil	1.35	0.31–0.29 (23°C) 0.30–0.28 (37°C)	0.15–0.19	✓	✓		✓		✓	✓	✓	

List of application notes used

- High-speed, high-resolution 3D imaging using resonant scanner and Denoise.ai
- The Advantages of Resonant Scanning with Ultra Short Laser Exposure Times in Live Imaging
- In vivo Confocal Imaging of Mouse Organs that Clearly Captures Fast Dynamics
- 3D Imaging of Intestinal Organoid



List of application notes



NIKON CORPORATION

Head office

1-5-20, Nishio, Shinagawa-ku, Tokyo, 140-8601, Japan
<https://www.healthcare.nikon.com/en/>

Manufacturer

471, Nagaodai-cho, Sakae-ku, Yokohama, Kanagawa 244-8533, Japan

ISO 14001 Certified
for NIKON CORPORATION

Specifications and equipment are subject to change without any notice or obligation on the part of the manufacturer.

July 2024 ©2023-2024 NIKON CORPORATION

N.B. Export of the products* in this brochure is controlled under the Japanese Foreign Exchange and Foreign Trade Law. Appropriate export procedures shall be required in case of export from Japan.

*Products: Hardware and its technical information (including software)

Nikon Instruments Inc.

1300 Walt Whitman Road, Melville, N.Y. 11747-3064, U.S.A.
 phone: +1-631-547-8500; +1-800-52-NIKON (within the U.S.A. only)
 fax: +1-631-547-0299

<https://www.microscope.healthcare.nikon.com/>

Nikon Europe B.V.

Stroombaan 14, 1181 VX Amstelveen, The Netherlands
 phone: +31-20-7099-000

https://www.microscope.healthcare.nikon.com/en_EU/

Nikon Precision (Shanghai) Co., Ltd.

CHINA phone: +86-21-6841-2050 fax: +86-21-6841-2060
 (Beijing branch) phone: +86-10-5831-2028 fax: +86-10-5831-2026
 (Guangzhou branch) phone: +86-2-3882-0551 fax: +86-2-3882-0580
<https://www.nikon-precision.com.cn/>

Nikon Canada Inc.

CANADA phone: +1-905-625-9910 fax: +1-905-602-9953

Nikon France, Succursale de Nikon Europe B.V.

FRANCE phone: +33-1-4516-4516

Nikon Deutschland, Zweigniederlassung der Nikon Europe B.V.

GERMANY phone: +49-211-9414-888

Nikon Italy, Branch of Nikon Europe B.V.

ITALY phone: +39-055-300-9601

Nikon Europe B.V., Amstelveen, Zweigniederlassung Schweiz (Egg/ZH)

SWITZERLAND phone: +41-43-277-2867

NIKON UK, Branch of Nikon Europe B.V.

UNITED KINGDOM phone: +44-208-247-1717

Nikon Österreich, Zweigniederlassung der Nikon Europe B.V.

AUSTRIA phone: +43-1-972-6111

Nikon Singapore Pte Ltd

SINGAPORE phone: +65-6559-3651 fax: +65-6559-3668

Nikon Australia Pty Ltd

AUSTRALIA phone: +61-2-8767-6900

Nikon Instruments Korea Co., Ltd.

KOREA phone: +82-2-6288-1900 fax: +82-2-555-4415

NIKON INDIA PVT. LTD.

AUSTRALIA phone: +91-124-4688-500



Code No. 2CE-MPHM-1 (2407) Am

WARNING

TO ENSURE CORRECT USAGE, READ THE CORRESPONDING MANUALS CAREFULLY BEFORE USING THE EQUIPMENT.